

Position Reconstruction of Rodent Movements based on Neural Spike Information of Place Cells

Günter Edlinger, Steve Schaffelhofer, Christoph Guger
g.tec medical engineering GmbH
Guger Technologies OEG
Graz, Austria
edlinger@gtec.at

J. Brotons-Mas, M. Sanchez-Vives
ICREA – Institut d'Investigacions Biomediques
August Pi i Sunyer
Barcelona, Spain

Abstract— Place cells are located in the hippocampus of the brain and play an important role for spatial navigation. In this study neural spike activity of freely moving rats along with the position of the rats was acquired. The study was performed to investigate if position reconstruction is possible if the rat is freely moving in open arenas with and without barriers of different rectangular sizes based on neural recordings. The neural spike activity was measured using tetrodes from 6 chronically implanted rats in CA1 and subiculum regions. The position of the rats was monitored using a video tracking system. In the encoding step spike activity features of the place cells and position information from the rats were used to train a computer system. In the decoding step, the position was estimated from the neural spiking activity. Different reconstruction methods were implemented: (i) Bayesian 1-step and (ii) Bayesian 2-step. The results show, that the reconstruction methods are able to determine the position of the rat in this 2-dimensional rectangular arenas from the cell activity. Better accuracies were achieved for CA1 regions because the firing fields show more localized spots. Furthermore, higher accuracies are achieved with a higher number of place cells and higher firing rates.

place cells, hippocampus, online position reconstruction, Bayesian 2-step method

I. INTRODUCTION

Cells with spatially modulated firing rates have been found in almost all areas of the hippocampus and in some surrounding areas. Such place cells (PC) encode an internal representation of an animal's position in its environment. The background firing rate of a PC is very low (0.1 Hz), but when an animal enters the receptive field of the neuron, its firing rate rapidly increases (up to ~ 5 - 15 Hz) [1, 2]. The location inducing the increased cell activity is called the firing field (FF) or place field (PF) of the particular neuron. However, also other sensory cues can influence place cell activity, but visual and motion-related cues are the most relevant [3].

Recently place cells were also used to reconstruct the foraging path of rats by investigating the firing field patterns [1, 5]. In an encoding stage the place cell spiking activity together with video tracking information is used to train a computer system on the data. In a decoding stage only the spike activity is used to reconstruct the path of the animal. In this study it was of interest to test if position reconstruction is also possible in open environments in contrast to the linear arenas used in previous studies [5]. The reconstruction was tested in square arenas with different side length (0,5m, 0,7m,

0,1m, 1m) and a square arena with a smaller square barrier inside (outer square: 1m x 1m, inner square: 0,6m x 0,6m). Two different algorithms based on Bayesian methods were implemented after a template matching approach was ruled out for reconstruction.

II. MATERIAL AND METHODS

A. Neural spike and video data acquisition

Action potential data were measured from 6 rats from CA1 or subiculum by the Instituto de Neurociencias de Alicante (Universidad Miguel Hernández-CSIC, UMH), Institute of Cognitive Neuroscience (University College London, UCL) and Center for Neural Systems, Memory and Aging (University of Arizona, UA). The rats were connected to the recording system via a head stage pre-amplifier. From 1 to 8 tetrodes were used for the recordings. Each data channel was amplified up to 40.000 times and a sampling frequency of 48.000 Hz was used. For tracking purpose small infrared light-emitting diodes were attached to the rat's head and a video camera system was mounted above the experimental arena. The number of recorded cells, the recording region and the arena sizes and shapes are shown in Table I. The sampling frequency for the position-tracking signal was 50 Hz (UMH, UA) and 46.875 Hz (UCL).

TABLE I. RECORDING INFORMATION OF THE 6 RATS.

Rat #	1	2	3	4	5	6
No. of cells	5	26	9	6	4	11
Hippocampal region	CA1	CA1	CA1	CA1	Subiculum	Subiculum
Test field side length [m]	0,7	0,8	0,5	1	1	1
Test field shape						

B. Dwell Time and Firing Rate

In the first step the recorded spike activity was used as input for the cluster cutting. The cluster cutting was performed manually to separate neuronal activity picked up by the electrodes into single cell activity. In a second step firing rate maps were created. Therefore the arena was divided into small subsets of space and class labels were assigned to these subsets

This work was supported by the EU-IST (FP6-027731) project Presencia, EU-IST (FP7-217148) project Synthetic Forager and EU-IST project Renachip. Special thanks to Paul Verschure and Matti Mintz for the discussions about place cells.

(pixels). Then the firing rate for each pixel was calculated by counting the number of spikes in each class. The test fields were divided into 64x64 classes, leading to an edge length for each class of 1.09 cm (rat 1), 1.25 cm (rat 2), 0.78 cm (rat 3) and 1.56 cm (rats 4, 5, 6). Figure 1A shows the movement trajectory of a rat obtained from the video tracking system for one full experiment with one rat.

The average firing rate of a cell i for each position \mathbf{x} is described by the firing rate distribution, i.e. a vector of firing rates per cell i : $f(\mathbf{x}) = (f_1(\mathbf{x}), f_2(\mathbf{x}), \dots, f_N(\mathbf{x}))$. In the training step, the firing rate maps $f(\mathbf{x})$ are calculated for the whole population of N recorded cells and for every single position \mathbf{x} . The average firing rates are calculated from the total number of spikes collected for a cell while the rat was at position $S(\mathbf{x})$, divided by the total amount of time the rat spent there $V(\mathbf{x})$ and Δt being the time interval of the position tracking system.

$$f(\mathbf{x}) = \frac{S(\mathbf{x})}{V(\mathbf{x}) \cdot \Delta t} \quad (1)$$

The firing rate distribution is independent of how often the rat populates a certain position, it rather describes the tendency of a cell to fire at a given position \mathbf{x} as shown in Fig. 1B for one neuron.

C. Position reconstruction algorithms

First firing rates of all cells within a sliding time window are compared with the firing rate vectors of each class that were setup in the training phase. Depending on the algorithm, a matching class number is returned which identifies the reconstructed position. To calculate the reconstruction error, the reconstructed position can be compared to the position known from the position tracking data. Then the time window is moved to the next reconstruction position.

For position reconstruction the two algorithms were implemented that were already tested in linear arenas: (i) Bayesian 1-step and (ii) Bayesian 2-step with continuity constraint. For the analysis all datasets were divided into three equally long parts. Training was performed on two thirds of the recordings; the other third was used for the testing and reconstruction.

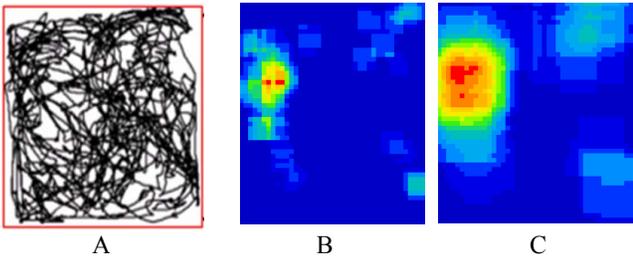


Figure 1. (A): example of movement trajectories of a rat recorded with the video tracking system seen from top. (B), (C): firing map of 1 neuron with different smoothing factors (5x5 and 10x10 Kernel).

1) Bayesian method (1-step)

When the rat moves to specific positions then certain spike trains are produced as response to the sensory stimuli. The probability of observing a number of spikes $n = (n_1, n_2, \dots, n_N)$ given the stimulus \mathbf{x} is $P(n|\mathbf{x})$. The probability of a stimulus (rat at position \mathbf{x}) to occur is defined as $P(\mathbf{x})$. The joint distribution

$$P(n, \mathbf{x}) = P(n|\mathbf{x}) \cdot P(\mathbf{x}) \quad (2)$$

measures the likelihood of observing both a stimulus \mathbf{x} and a spike train n .

For the reconstruction, the observation of a number of spikes n found in the measured data occur with a probability of $P(n)$. The likelihood of observing both stimulus and spikes is finally given by [5]

$$P(\mathbf{x} | n) = C(\tau, n) \cdot P(\mathbf{x}) \cdot \left(\prod_{i=1}^N f_i(\mathbf{x})^{n_i} \right) \cdot \exp\left(-\tau \cdot \sum_{i=1}^N f_i(\mathbf{x})\right) \quad (3)$$

where $C(\tau, n)$ is a normalization factor and τ is the length of the time window. In this study the normalization factor $C(\tau, n)$ was set to 1.

The most probable position will be regarded as the animal's position:

$$\hat{\mathbf{x}}_{Bayes, 1-step} = \arg \max_{\mathbf{x}} P(\mathbf{x} | n) \quad (4)$$

2) Bayesian method with continuity constraint (2-step)

A continuity constraint can improve the accuracy of reconstruction. Sudden jumps in the reconstruction are observed by the low instantaneous firing rates of the recorded place cells. If not enough cells are firing then there is a lack of information. However the firing information is needed for the position reconstruction. The continuity constraint incorporates the reconstructed position from the preceding time step as well as the current speed information.

Based on the 1-step Bayesian method the reconstructed position at the preceding time step $t-1$ is also used to calculate the reconstructed position at time t :

$$P(\mathbf{x}_t | n_t, \mathbf{x}_{t-1}) = k \cdot P(\mathbf{x}_t | n_t) \cdot P(\mathbf{x}_{t-1} | \mathbf{x}_t) \quad (5)$$

For more detail on this approach see [5].

III. RESULTS

Fig. 2 shows the reconstruction results computed with the Bayesian 2-Step algorithm of rat 3. Data from seconds 160 to 480 were used to train the algorithm and the interval 0 to 160

was used to test the method. The figure shows that the reconstructed path follows the real path very well in many data points of the recording. The mean reconstruction error of rat 3 was 9.5 cm. Fig. 2 shows clearly erratic jumps of the reconstructed path in both x- and y-coordinates.

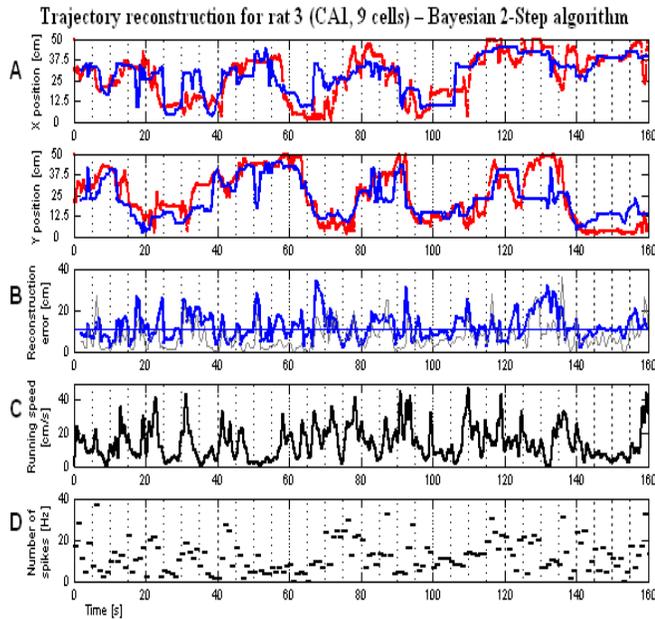


Figure 2. Rat's 3 reconstruction trajectory based on place cell firing. (A): Real (red) and reconstructed (blue) x- and y-positions. (B): Reconstruction error (thick line), reconstruction error weighted by the instant place cell firing rate (thin line) and median error of the whole recording (horizontal line). (C): Running speed. (D): Firing rate of all neurons.

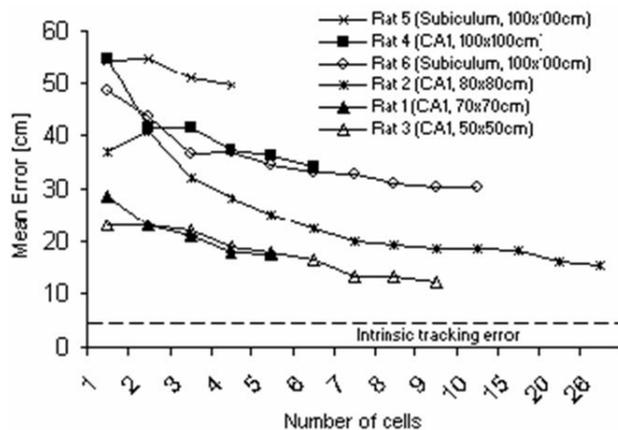


Figure 3. Mean reconstruction error for all rats as a function of number of cells.

If the reconstruction is only performed for time intervals where a minimum number of spikes (4 spikes in reconstruction window) are present, the accuracy is increased as shown by the thin gray thin line in Figure 2B. In this case

the mean error is 8 cm. Interesting is the erratic jump around second 51 where the running speed (Fig. 2C) was almost 0.

In Figure 3 the mean reconstruction error based on the Bayesian 2-step algorithm for all six rats are given as a function of the number of recorded cells. The horizontal dotted line at 5cm in the graph displays the intrinsic tracking error, which is defined as the average uncertainty in position tracking, due to the size of the diode (LED) arrays, the distance of the diodes above the rat's head and variations in posture [6]. Here it can be seen that for rat 3 the mean reconstruction error drops from 22 cm for 1 unit to about 10 cm for 9 units.

IV. DISCUSSION

The reconstruction algorithms for neural spike trains were implemented and applied to place cell activity from the hippocampus. The goal was to reconstruct the position of the rats just by investigating the spike activity as accurate as possible and therefore to minimize the reconstruction error. The best performance was found for the Bayesian 2-Step algorithm. The reason is that the algorithm considers also the previous position of the rat and does not allow large jumps from one reconstruction point to the next one. The Bayesian 1-Step algorithm performed less accurate but it is very interesting to see the results because it is only based on the current reconstruction window.

Despite subicular units tend to be more stable than CA1 cells, the results show that position reconstruction is more accurate with CA1 units. However, it is also interesting that subicular units contain enough information for the reconstruction. This can also be seen from the place field density plots, where the place fields of subicular units are much more blurred than place fields of CA1 units. Distinct place fields in combination with a high number of cells guarantee good reconstruction results. But it must be noted that only 2 data-sets from the subiculum and 4 data-sets from CA1 regions were investigated. The investigation of more data-sets is necessary to prove that.

Different accuracies in position reconstruction have been reported. Zhang reached classification accuracies of about 4-10 cm with about 30 cells in 2 rats [5]. Wilson [7] reported an error of 33 cm with 10 cells. Barbieri [8] achieved an error of 5.1 cm in a 70 cm circular environment with 44 neurons. Theoretically the reconstruction error is inversely proportional to the square root of the number of cells. This has been confirmed by the data analysis of all 6 rats as well as in other publications [5, 7]. Interesting is that for 3 rats (1, 2, 3) the reconstruction error was below 20 cm with less than 10 place cells. This study aimed to investigate bigger square shaped arenas. Results show that with even a few cells the reconstruction can already be performed. But important is to consider the different arena sizes.

Erratic jumps occur also when the animal stops running. This has two reasons: (i) the firing rate is modulated by speed and therefore lower and (ii) the animal received food rewards to move around in the arena and therefore the eating can

produce artifacts in the neural data recorded. Zhang also suggests biological reasons for the erratic jumps like the rat is looking around or planning the next move [5]. Another explanation of the erratic jumps in the reconstruction could be the forward and reverse replay of cell firing during immobility of the rat [9]. As already mentioned in the results section, minimizing the reconstruction to points where the firing rate is above a certain threshold only punctually yields to better reconstruction results. It does not reduce the overall error rate, because the positions that have not been reconstructed have to be interpolated. This estimation leads to reconstruction errors for the interpolated positions again, and increases the error rate.

It must be noted that the behavior of the rodents in square shaped arenas might be also different from the linear and circular environments studied earlier as the rodents like to create a home base in corners. Eilam investigated home base behavior and found that excursions starting from a home base needs longer times than coming back [10].

Next steps in this research will include the investigation of grid cells, the consideration of the accuracy with respect to the taken path with respect to the home base as well as the realization of real-time reconstruction hardware and software setup.

REFERENCES

- [1] O'Mara S. M., "Spatially selective firing properties of hippocampal formation neurons in rodents and primates," *Progress in Neurobiology* Vol. 45, 253-274, 1995.
- [2] Anderson M. I. and O'Mara S. M., "Analysis of recordings of single-unit firing and population activity in the dorsal subiculum of unrestrained, freely moving rats," *The Journal of Neurophysiology* Vol. 90 No. 2, 655-665, 2003.
- [3] Poucet B., Lenck-Santini P. P., Paz-Villagrán V. and Save E., "Place cells, neocortex and spatial navigation: a short review," *J Physiol Paris* Vol. 97 Issues 4-6, 537-546, 2003.
- [4] Brown E. N., Frank L. M., Tang D., Quirk M. C. and Wilson M. A., "A Statistical Paradigm for Neural Spike Train Decoding Applied to Position Prediction from Ensemble Firing Patterns of Rat Hippocampal Place Cells," *The Journal of Neuroscience* Vol. 18(18), 7411-7425, 1998.
- [5] Zhang K., Ginzburg I., McNaughton B. L. and Sejnowski T. J., "Interpreting Neuronal Population Activity by Reconstruction: Unified Framework With Application to Hippocampal Place Cells," *The Journal of Neurophysiology* Vol. 79, No. 2, 1017-1044, 1998.
- [6] Skaggs W. E., McNaughton B. L., Wilson M. A. and Barnes C. A., "Theta Phase Precession in Hippocampal Neuronal Populations and the Compression of Temporal Sequences," *Hippocampus* Vol. 6(2), 149-172, 1996.
- [7] Wilson M. A. and McNaughton B. L., "Dynamics of the Hippocampal Ensemble Code for Space," *Science* Vol. 261 (5124), 1055-1058, 1993.
- [8] Barbieri R., M.A. Wilson, L.M. Frank, and Brown E.N., "An analysis of hippocampal spatio-temporal representations using a Bayesian algorithm for neural spike train decoding," *IEEE Trans. Neural Systems and Rehab. Engng.* Vol.13(2), 131-136, 2005.
- [9] Diba K. and Buzsaki G., "Forward and reverse hippocampal place-cell sequences during ripples," *Nat. Neuroscience* 10:1241-2, 2007.
- [10] Eilam D. and Golani I., "Home base behavior of rats (*Rattus norvegicus*) exploring a novel environment," *Behav Brain Res.*, 34(3):199-211, 1989.